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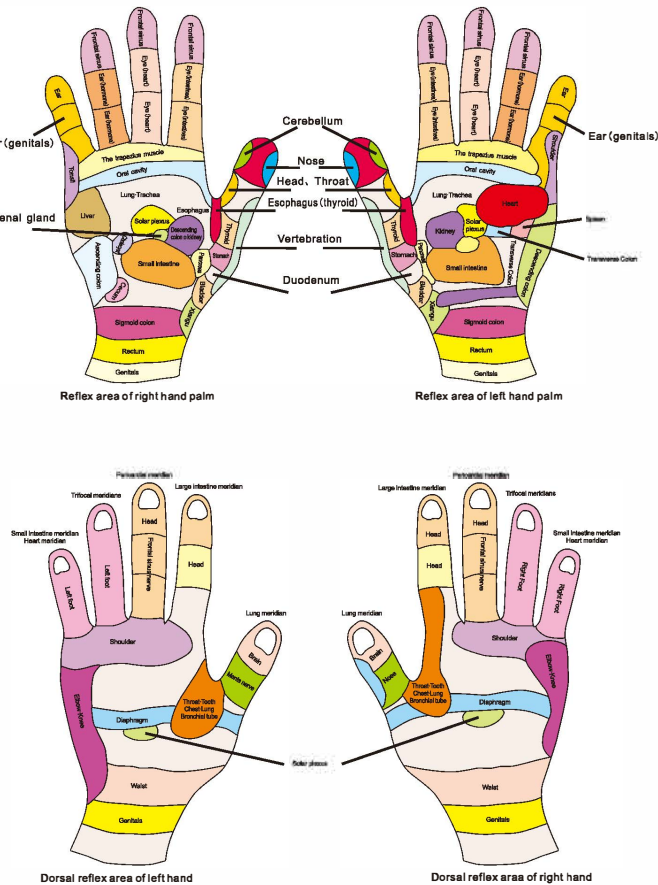
Thank you for choosing this product. Before you use it, please read the instructions in order to master the correct use of this product. When referring, please pay special attention to the matters needing attention, and keep the instructions properly for reference.

Note: We reserve the right to make design changes and interpretations of our products. We will not inform you if there are any changes. Product color is for reference only, please refer to the actual object.

Content

1、Product description	1
2、Operation instructions	2
3、Usage method	3
4、Attentions	4
5、Maintenance and storage	5
6、Product parameters	5
7、Troubleshooting	6
8、Palm acupuncture point	7


If there is a malfunction in the operation of this product, please check it according to the following methods. If the malfunction can not be eliminated or the product is still unable to operate, please turn off the power and contact the product service center for maintenance.

Basic Faults	Possible causes	Troubleshooting Measures
The product can't start working	The remaining power may be weak.	Charging
The product cannot be recharged	Is the power adapter connected to the data line correctly?	Re-plug until LCD power display flickers
	Is the data line and USB socket intact?	Check the socket. If there is any damage, please contact the manufacturer.
	Are the specifications of adapters and data lines in line with operational requirements?	View specifications and proofread the technical parameters provided in the instructions
The product can be turned on, but can't work	The remaining power may be weak.	Charging

Do not disassemble or attempt to repair by yourself

1. Use towel to remove dust and other dirt, and place it in a dry and ventilated place.
2. Do not place in high temperature, humidityMaintenance and storage and direct sunlight.
3. When the storage period exceeds 6 months, it should be checked out and power-on.

Product name : Hand massager
Rated power:8W
Rated voltage: 5V~1000mA
Protection time: 5minutes、10minutes、15minutes
Battery capacity:2200mAh
Power type:lithium battery
Massage mode:Air bag
Size: 308X158X116mm
Product weight:1kg
Executive standard:GB4706.1-2005 GB4706.10 2008

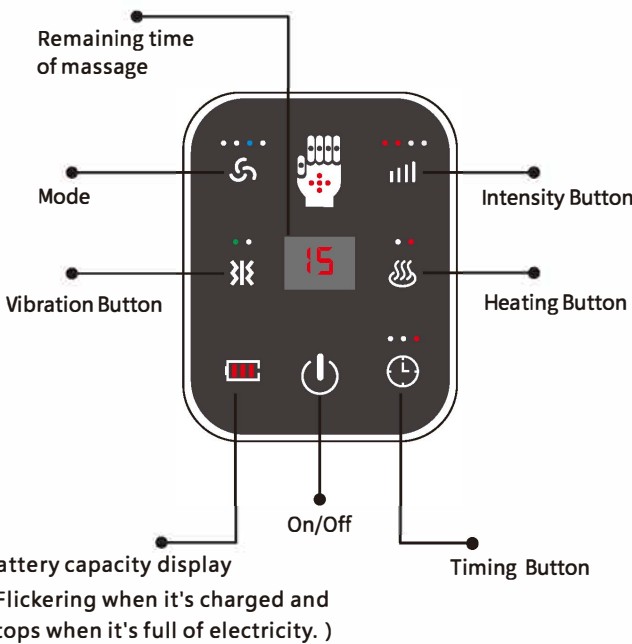1. Put your hand inside the machine.
2. Put the palm in the center of the massage pad.
3. Press the On/Off button (⏻).
4. Select the massage mode (🌀): yellow light (vitality mode), green light (finger+ palm massage), blue light (wrist massage), and red light (deep massage).
5. Select work intensity, (||||) set intensity, a total of 4 work intensities.
6. If you need to turn on the heating function, please press (🔥) to turn on the heating mode. The green light is low temperature, the red light is high temperature.
7. If you need to turn on the vibration, please select the button (🌀), there are two different vibration modes.
8. Timing setting: Press the button (⌚) to select 3 types of timing: 5, 10, and 15 minutes. (You can check the remaining working time through the LCD display)
9. Please press the button (⏻) to turn off the machine.

Tips :

- ★Please don't use products with heating function for a long time to prevent low temperature burns.
- ★When the product is charged, the power sign flickers until it stops when it is full.
- ★Massage with plastic gloves after applying hand cream can help to absorb hand cream and repair skin.

1. Use the type-C plug connection. (Note the shape of the insertion port)
2. Please make sure that there are no jewelry (watch ring Bracelet ring, etc.) in your hands before use so as to avoid the damage of jewelry caused by sharp objects piercing the airbag or squeezing the airbag.
3. It is forbidden to use for a long time. Normally, massage time should not exceed 30 minutes.
4. Please do not use this product while driving.
5. Please do not use it after getting drunk or taking sleeping pills.
6. Please use our charger (optional) or DC5V 1 ~ 4A charger.
7. This product is designed for household electronic health care products. It is not suitable for diagnosis treatment and other professional or commercial purposes.
8. Please keep this product away from water, fire and corrosive environment.
9. Please do not put sharp objects into this product.
10. Please do not use when hands are wet.
11. If the power cord is damaged stop using it immediately.
12. If you feel unwell, please stop using it in time.
13. Please do not wrap products with towels, quilts, clothes, etc. when you use it.

Please follow the doctor's advice for the following patients :
The patients who have had hand surgery, swelling and inflammation, congestion, trauma and other hand diseases;
Blood diseases with bleeding tendency; osteoporosis and soft tissue diseases; infants and pregnant women;